Epstein-Barr Virus (VCA) IgM

ELISA

CE

Only for in-vitro diagnostic use

English	2
Deutsch	7
Français	13
Italiano	19
Español	24
Português.....	29
Bibliography / Literatur / Bibliographie / Bibliografia / Bibliografía/ Bibliografia	34
Abbreviations / Abkürzungen / Abréviations / Abbreviazioni / Abreviaciones / Abreviaturas	34
Symbols Key / Symbolschlüssel / Explication des Symboles / Legenda / Símbolos / Tabela de símbolos	35
Summary of Test Procedure / Kurzanleitung Testdurchführung / Résumé de la procédure de test / Schema della procedura / Resumen de la técnica / Resumo do Procedimento de Teste	36

Product Number: EBVM0150BA (96 Determinations)

ENGLISH

1. INTRODUCTION

Epstein-Barr Virus (EBV) is a member of the herpes virus family (Gamma subgroup, DNA virus of 120-200 nm) and one of the most common human viruses. The virus occurs worldwide, and most people become infected with EBV sometime during their lives. Transmission of the virus is almost impossible to prevent since many healthy people can carry and spread the virus intermittently for life. Infants become susceptible to EBV as soon as maternal antibody protection disappears. Infection of children usually causes no symptoms. Infection during adolescence or young adulthood causes infectious mononucleosis 35% to 50% of the time.

Infectious mononucleosis is almost never fatal. There are no known associations between active EBV infection and problems during pregnancy, such as miscarriages or birth defects. Although the symptoms of infectious mononucleosis usually resolve in 1 or 2 months, EBV remains dormant or latent in a few cells in the throat and blood for the rest of the person's life. Periodically, the virus can reactivate and is commonly found in the saliva of infected persons. This reactivation usually occurs without symptoms of illness.

EBV also establishes a lifelong dormant infection in some cells of the body's immune system. A late event in a very few carriers of this virus is the emergence of Burkitt's lymphoma and nasopharyngeal carcinoma, but EBV is probably not the sole cause of these malignancies.

Species	Disease	Symptoms (e.g.)	Transmission route
Epstein-Barr Virus	Infectious mononucleosis	Fever, sore throat, swollen lymph glands	Person to Person transmission by saliva

Infection or presence of pathogen may be identified by:

- PCR
- Serology: "mono spot" test, Detection of antibodies by ELISA

2. INTENDED USE

The Epstein-Barr Virus (VCA) IgM ELISA is intended for the qualitative determination of IgM class antibodies against Epstein-Barr Virus **viral capsid antigen** (VCA) in human serum or plasma (citrate, heparin).

3. PRINCIPLE OF THE ASSAY

The qualitative immunoenzymatic determination of specific antibodies is based on the ELISA (Enzyme-linked Immunosorbent Assay) technique.

Microtiterplates are coated with specific antigens to bind corresponding antibodies of the sample. After washing the wells to remove all unbound sample material a horseradish peroxidase (HRP) labelled conjugate is added. This conjugate binds to the captured antibodies. In a second washing step unbound conjugate is removed. The immune complex formed by the bound conjugate is visualized by adding Tetramethylbenzidine (TMB) substrate which gives a blue reaction product.

The intensity of this product is proportional to the amount of specific antibodies in the sample. Sulphuric acid is added to stop the reaction. This produces a yellow endpoint colour. Absorbance at 450/620 nm is read using an ELISA Microtiterplate reader.

4. MATERIALS

4.1. Reagents supplied

- **Microtiterplate:** 12 break-apart 8-well snap-off strips coated with Epstein-Barr Virus (VCA) synthetic p18 peptide; in resealable aluminium foil.
- **IgM Sample Dilution Buffer:** 1 bottle containing 100 mL of phosphate buffer (10 mM) for sample dilution; pH 7.2 ± 0.2; anti-human IgG (RF Absorbent); coloured green; ready to use; white cap; ≤ 0.0015% (v/v) CMIT/ MIT (3:1).
- **Stop Solution:** 1 bottle containing 15 mL sulphuric acid, 0.2 mol/L; ready to use; red cap.
- **Washing Buffer (20x conc.):** 1 bottle containing 50 mL of a 20-fold concentrated phosphate buffer (0.2 M), pH 7.2 ± 0.2, for washing the wells; white cap.
- **Conjugate:** 1 bottle containing 20 mL of peroxidase labelled antibody to human IgM in phosphate buffer (10 mM); coloured red; ready to use; black cap.
- **TMB Substrate Solution:** 1 bottle containing 15 mL 3,3',5,5'-tetramethylbenzidine (TMB), < 0.1%; ready to use; yellow cap.
- **Positive Control:** 1 vial containing 2 mL control; coloured yellow; ready to use; red cap; ≤ 0.02% (v/v) MIT.
- **Cut-off Control:** 1 vial containing 3 mL control; coloured yellow; ready to use; green cap; ≤ 0.02% (v/v) MIT.
- **Negative Control:** 1 vial containing 2 mL control; coloured yellow; ready to use; blue cap; ≤ 0.0015% (v/v) CMIT/ MIT (3:1).

For hazard and precautionary statements see 12.1

For potential hazardous substances please check the safety data sheet.

4.2. Materials supplied

- 1 Cover foil
- 1 Instruction for use (IFU)
- 1 Plate layout

4.3. Materials and Equipment needed

- ELISA Microtiterplate reader, equipped for the measurement of absorbance at 450/620 nm
- Incubator 37 °C
- Manual or automatic equipment for rinsing Microtiterplates
- Pipettes to deliver volumes between 10 and 1000 µL
- Vortex tube mixer
- Distilled water
- Disposable tubes

5. STABILITY AND STORAGE

Store the kit at 2...8 °C. The opened reagents are stable up to the expiry date stated on the label when stored at 2...8 °C.

6. REAGENT PREPARATION

It is very important to bring all reagents and samples to room temperature (20...25 °C) and mix them before starting the test run!

6.1. Microtiterplate

The break-apart snap-off strips are coated with Epstein-Barr Virus (VCA) synthetic p18 peptide. Immediately after removal of the strips, the remaining strips should be resealed in the aluminium foil along with the desiccant supplied and stored at 2...8 °C.

6.2. Washing Buffer (20x conc.)

Dilute Washing Buffer 1 + 19; e. g. 10 mL Washing Buffer + 190 mL distilled water. The diluted buffer is stable for 5 days at room temperature (20...25 °C). In case crystals appear in the concentrate, warm up the solution to 37 °C e.g. in a water bath. Mix well before dilution.

6.3. TMB Substrate Solution

The reagent is ready to use and has to be stored at 2...8 °C, away from the light. The solution should be colourless or could have a slight blue tinge. If the substrate turns into blue, it may have become contaminated and should be thrown away.

7. SAMPLE COLLECTION AND PREPARATION

Use human serum or plasma (citrate, heparin) samples with this assay. If the assay is performed within 5 days after sample collection, the samples should be kept at 2...8 °C; otherwise they should be aliquoted and stored deep-frozen (-70...-20 °C). If samples are stored frozen, mix thawed samples well before testing. Avoid repeated freezing and thawing. Heat inactivation of samples is not recommended.

7.1. Sample Dilution

Before assaying, all samples should be diluted 1+100 with IgM Sample Dilution Buffer. Dispense 10 µL sample and 1 mL IgM Sample Dilution Buffer into tubes to obtain a 1+100 dilution and thoroughly mix with a Vortex.

8. ASSAY PROCEDURE

Please read the instruction for use carefully **before** performing the assay. Result reliability depends on strict adherence to the instruction for use as described. The following test procedure is only validated for manual procedure. If performing the test on ELISA automatic systems we recommend increasing the washing steps from three up to five and the volume of Washing Buffer from 300 µL to 350 µL to avoid washing effects. Pay attention to chapter 12. Prior to commencing the assay, the distribution and identification plan for all samples and standards/controls (duplicates recommended) should be carefully established on the plate layout supplied in the kit. Select the required number of microtiter strips or wells and insert them into the holder.

Perform all assay steps in the order given and without any delays.

A clean, disposable tip should be used for dispensing each standard/control and sample.

Adjust the incubator to 37 ± 1 °C.

1. Dispense 100 µL standards/controls and diluted samples into their respective wells. Leave well A1 for the Substrate Blank.
2. Cover wells with the foil supplied in the kit.
3. **Incubate for 1 hour \pm 5 min at 37 ± 1 °C.**
4. When incubation has been completed, remove the foil, aspirate the content of the wells and wash each well three times with 300 µL of Washing Buffer. Avoid overflows from the reaction wells. The interval between washing and aspiration should be > 5 sec. At the end carefully remove remaining fluid by tapping strips on tissue paper prior to the next step!
Note: Washing is important! Insufficient washing results in poor precision and false results.
5. Dispense 100 µL Conjugate into all wells except for the Substrate Blank well A1.
6. **Incubate for 30 min at room temperature (20...25 °C).** Do not expose to direct sunlight.
7. Repeat step 4.
8. Dispense 100 µL TMB Substrate Solution into all wells.
9. **Incubate for exactly 15 min at room temperature (20...25 °C) in the dark.** A blue colour occurs due to an enzymatic reaction.
10. Dispense 100 µL Stop Solution into all wells in the same order and at the same rate as for the TMB Substrate Solution, thereby a colour change from blue to yellow occurs.
11. Measure the absorbance at 450/620 nm within 30 min after addition of the Stop Solution.

8.1. Measurement

Adjust the ELISA Microtiterplate reader **to zero** using the **Substrate Blank**.

If - due to technical reasons - the ELISA Microtiterplate reader cannot be adjusted to zero using the Substrate Blank, subtract its absorbance value from all other absorbance values measured in order to obtain reliable results!

Measure the absorbance of all wells at **450 nm** and record the absorbance values for each standard/control and sample in the plate layout.

Bichromatic measurement using a reference wavelength of 620 nm is recommended.

Where applicable calculate the mean absorbance values of all duplicates.

9. RESULTS

9.1. Run Validation Criteria

In order for an assay run to be considered valid, these Instructions for Use have to be strictly followed and the following criteria must be met:

- **Substrate Blank:** Absorbance value **< 0.100**
- **Negative Control:** Absorbance value **< 0.200 and $< \text{Cut-off}$**
- **Cut-off Control:** Absorbance value **$0.150 - 1.300$**
- **Positive Control:** Absorbance value **$> \text{Cut-off}$**

If these criteria are not met, the test is not valid and must be repeated.

9.2. Calculation of Results

The Cut-off is the mean absorbance value of the Cut-off Control determinations.

Example: Absorbance value Cut-off Control 0.44 + absorbance value Cut-off control 0.42 = 0.86 / 2 = 0.43
Cut-off = 0.43

9.2.1. Results in Units [U]

$\frac{\text{Sample (mean) absorbance value} \times 10}{\text{Cut-off}} = [\text{Bioactive Units} = \text{U}]$

Example: $\frac{1.591 \times 10}{0.43} = 37 \text{ U (Units)}$

9.3. Interpretation of Results

Cut-off	10 U	-
Positive	> 11 U	Antibodies against the pathogen are present. There has been a contact with the antigen (pathogen resp. vaccine).
Equivocal	9 – 11 U	Antibodies against the pathogen could not be detected clearly. It is recommended to repeat the test with a fresh sample in 2 to 4 weeks. If the result is equivocal again the sample is judged as negative .
Negative	< 9 U	The sample contains no antibodies against the pathogen. A previous contact with the antigen (pathogen resp. vaccine) is unlikely.
Diagnosis of an infectious disease should not be established on the basis of a single test result. A precise diagnosis should take into consideration clinical history, symptomatology as well as serological data. In immunocompromised patients and newborns serological data only have restricted value.		

9.3.1. Antibody Isotypes and State of Infection

Interpretation of results depends on the specific clinical application of the test: any laboratory should establish its own clinically relevant ranges for the population taken into consideration. Prevalence may very depending on geographical location, age, socioeconomic status, type of test employed, specimen collection and handling procedures, clinical and epidemiological history of individual patients.

Antibody profile of EBV infections			Stage of EBV infection
VCA IgM	VCA IgG	EBNA IgG	
-	-	-	EBV negative ¹
+	-	-	Primary EBV infection (early phase) ²
+	+	-	Primary EBV infection (acute phase)
-	+	-	Uncertain result ³
-	+	+	Past EBV infection
-	-	+	Uncertain result ⁴
+	+	+	Uncertain result ⁵

¹ on suspicion of virus exposition, examination of a second sample ca. 7 days later

² clarification: cross reactive IgM in combination with CMV primary infection

³ differentiation: primary infection with lack of IgM or past infection with negative EBNA1 necessary

⁴ very rare constellation of a past infection; clarification of unspecific EBNA 1 necessary

⁵ differentiation: shortly past primary infection, reactivation, cross reactive IgM of CMV primary infection or polyclonal IgM stimulation necessary

10. SPECIFIC PERFORMANCE CHARACTERISTICS

The results refer to the groups of samples investigated; these are not guaranteed specifications.

10.1. Precision

Intraassay	n	Mean (E)	CV (%)
#1	24	0.502	2.28
#2	24	1.542	1.25
#3	24	0.836	3.18
Interassay	n	Mean (U)	CV (%)
#1	12	30.07	6.89
#2	12	13.73	9.40
#3	12	2.26	9.30

10.2. Diagnostic Specificity

The diagnostic specificity is defined as the probability of the assay of scoring negative in the absence of the specific analyte. It is 97.62% (95% confidence interval: 94.02% - 99.35%).

10.3. Diagnostic Sensitivity

The diagnostic sensitivity is defined as the probability of the assay of scoring positive in the presence of the specific analyte. It is 94.59% (95% confidence interval: 81.81% - 99.34%).

10.4. Interferences

Interferences with hemolytic, lipemic or icteric samples are not observed up to a concentration of 10 mg/mL hemoglobin, 5 mg/mL triglycerides and 0.5 mg/mL bilirubin.

10.5. Cross Reactivity

Investigation of a sample panel with antibody activities to potentially cross-reacting parameters did not reveal evidence of false-positive results due to cross-reactions.

11. LIMITATIONS OF THE PROCEDURE

Bacterial contamination or repeated freeze-thaw cycles of the sample may affect the absorbance values.

12. PRECAUTIONS AND WARNINGS

- The test procedure, the information, the precautions and warnings in the instructions for use have to be strictly followed. The use of the testkits with analyzers and similar equipment has to be validated. Any change in design, composition and test procedure as well as for any use in combination with other products not approved by the manufacturer is not authorized; the user himself is responsible for such changes. The manufacturer is not liable for false results and incidents for these reasons. The manufacturer is not liable for any results by visual analysis of the patient samples.
- Only for in-vitro diagnostic use.
- All materials of human or animal origin should be regarded and handled as potentially infectious.
- All components of human origin used for the production of these reagents have been tested for anti-HIV antibodies, anti-HCV antibodies and HBsAg and have been found to be non-reactive.
- Do not interchange reagents or Microtiterplates of different production lots.
- No reagents of other manufacturers should be used along with reagents of this test kit.
- Do not use reagents after expiry date stated on the label.
- Use only clean pipette tips, dispensers, and lab ware.
- Do not interchange screw caps of reagent vials to avoid cross-contamination.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- After first opening and subsequent storage check conjugate and standard/control vials for microbial contamination prior to further use.
- To avoid cross-contamination and falsely elevated results pipette patient samples and dispense reagents without splashing accurately into the wells.
- The ELISA is only designed for qualified personnel following the standards of good laboratory practice (GLP).
- For further internal quality control each laboratory should additionally use known samples.

12.1. Safety note for reagents containing hazardous substances

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Therefore, the following hazard and precautionary statements apply.

	Warning	H317	May cause an allergic skin reaction.
		P261	Avoid breathing spray.
		P280	Wear protective gloves/ protective clothing.
		P302+P352	IF ON SKIN: Wash with plenty of soap and water.
		P333+P313	If skin irritation or rash occurs: Get medical advice/ attention.
		P362+P364	Take off contaminated and Wash it before reuse.

Further information can be found in the safety data sheet

12.2. Disposal Considerations

Residues of chemicals and preparations are generally considered as hazardous waste. The disposal of this kind of waste is regulated through national and regional laws and regulations. Contact your local authorities or waste management companies which will give advice on how to dispose hazardous waste.

13. ORDERING INFORMATION

Prod. No.: EBVM0150BA Epstein-Barr Virus (VCA) IgM ELISA (96 Determinations)

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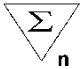
- Bergman, Michael M.; Gleckman, Richard A. (1987): Heterophil-negative infectious mononucleosis-like syndrome. In *Postgraduate medicine* 81 (1), 313-20, 325-6.
- Buchwald, D.; Komaroff, A. L. (1991): Review of laboratory findings for patients with chronic fatigue syndrome. In *Reviews of infectious diseases* 13 Suppl 1, S12-8.
- De-The, G. (1983): Epidemiology of EBV and associated diseases in man. In Heinz Fraenkel-Conrat, Robert R. Wagner, Bernard Roizman (Eds.): *Herpesviruses*. N. Boston: Springer US, pp. 25–103.
- Dölken, G.; Weitzmann, U.; Boldt, C.; Bitzer, M.; Brugger, W.; Lohr, G. W. (1984): Enzyme-linked immunosorbent assay for IgG antibodies to Epstein-Barr virus-associated early antigens and viral capsid antigen. In *Journal of immunological methods* 67 (2), pp. 225–233.
- Färber, Inge; Wutzler, Peter; Wohlrabe, Peter; Wolf, Hans; Hinderer, Walter; Sonneborn, Hans-Hermann (1993): Serological diagnosis of infectious mononucleosis using three anti-Epstein-Barr virus recombinant ELISAs. In *Journal of virological methods* 42 (2-3), pp. 301–307.
- Gorgievski-Hrisoho, Meri; Hinderer, Walter; Nebel-Schickel, Heike; Horn, Jürgen; Vornhagen, Rolf; Sonneborn, Hans-Hermann et al. (1990): Serodiagnosis of infectious mononucleosis by using recombinant Epstein-Barr virus antigens and enzyme-linked immunosorbent assay technology. In *Journal of Clinical Microbiology* 28 (10), pp. 2305–2311.
- Halprin, Jessica; Scott, Alan L.; Jacobson, Lisa; Levine, Paul H.; Ho, J. H.C.; Niederman, James C. et al. (1986): Enzyme-linked immunosorbent assay of antibodies to Epstein-Barr virus nuclear and early antigens in patients with infectious mononucleosis and nasopharyngeal carcinoma. In *Annals of internal medicine* 104 (3), pp. 331–337.
- Heath, Clark W., JR; Brodsky, Allan L.; Potolsky, Abraham I. (1972): Infectious mononucleosis in a general population. In *American Journal of Epidemiology* 95 (1), pp. 46–52.
- Henle, Werner; Henle, Gertrude E.; Horwitz, Charles A. (1974): Epstein-Barr virus specific diagnostic tests in infectious mononucleosis. In *Human pathology* 5 (5), pp. 551–565.
- Hudnall, S. David; Stanberry, Lawrence R. (2006): Human Herpesvirus Infections. In Richard L. Guerrant, David H. Walker, Peter F. Weller (Eds.): *Tropical infectious diseases. Principles, pathogens & practice*. 2nd ed. Philadelphia: Churchill Livingstone, pp. 590–620.
- Lamy, M. E.; Favart, A. M.; Cornu, C.; Mendez, M.; Segas, M.; Burtonboy, G. (1982): Study of Epstein Barr virus (EBV) antibodies: IgG and IgM anti-VCA, IgG anti-EA and Ig anti-EBNA obtained with an original microtiter technique: --serological criterions of primary and recurrent EBV infections and follow-up of infectious mononucleosis--seroepidemiology of EBV in Belgium based on 5178 sera from patients. In *Acta clinica Belgica* 37 (5), pp. 281–298.
- Lennette, E. (1985): Epstein-Barr Virus. In : *Manual of clinical microbiology*. Edited by Edwin H. Lennette. Washington: American Society for microbiology, pp. 728–732.
- Luka, Janos; Chase, Robert C.; Pearson, Gary R. (1984): A sensitive enzyme-linked immunosorbent assay (ELISA) against the major EBV-associated antigens. I. Correlation between ELISA and immunofluorescence titers using purified antigens. In *Journal of immunological methods* 67 (1), pp. 145–156.
- Ory, Fernando de; Antonaya, José; Fernández, María Visitación; Echevarría, José Manuel (1993): Application of low-avidity immunoglobulin G studies to diagnosis of Epstein-Barr virus infectious mononucleosis. In *Journal of Clinical Microbiology* 31 (6), pp. 1669–1671.
- Pearson, G. P. (2012): Infectious Mononucleosis: the humoral response. In H.-Harald Sedlacek (Ed.): *Infectious mononucleosis*. [Place of publication not identified]: Springer.
- Pochedly, Carl (1987): Laboratory testing for infectious mononucleosis. Cautions to observe in interpreting results. In *Postgraduate medicine* 81 (1), 335-9, 342.

ABBREVIATIONS / ABKÜRZUNGEN / ABRÉVIATIONS / ABBREVIAZIONI / ABREVIACIONES / ABREVIATURAS

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Epstein-Barr Virus (VCA) IgM ELISA

Test Preparation

Prepare reagents and samples as described.
 Establish the distribution and identification plan for all samples and standards/controls on the plate layout supplied in the kit.
 Select the required number of microtiter strips or wells and insert them into the holder.

Assay Procedure

	Substrate Blank (A1)	Negative Control	Cut-off Control	Positive Control	Sample (diluted 1+100)
Negative Control	-	100 µL	-	-	-
Cut-off Control	-	-	100 µL	-	-
Positive Control	-	-	-	100 µL	-
Sample (diluted 1+100)	-	-	-	-	100 µL
Cover wells with foil supplied in the kit Incubate for 1 h at 37 ± 1 °C Wash each well three times with 300 µL of Washing Buffer					
Conjugate	-	100 µL	100 µL	100 µL	100 µL
Incubate for 30 min at room temperature (20...25 °C) Do not expose to direct sunlight Wash each well three times with 300 µL of Washing Buffer					
TMB Substrate solution	100 µL	100 µL	100 µL	100 µL	100 µL
Incubate for exactly 15 min at room temperature (20...25 °C) in the dark					
Stop Solution	100 µL	100 µL	100 µL	100 µL	100 µL
Photometric measurement at 450 nm (reference wavelength: 620 nm)					

Bioactiva Diagnostica GmbH

Louisenstrasse 137
 61348 Bad Homburg, Germany

Tel: +49 (0) 6172 17102-0 Fax: +49 (0) 6172 17102-29

Email: bioactiva@bioactiva.de

Internet: www.bioactiva.de

EBVM0150BA_2020-08-09_Ka-ab Lot 184